

OriGene Technologies, Inc.


9620 Medical Center Drive, Ste 200 Rockville, MD 20850, US Phone: +1-888-267-4436 https://www.origene.com techsupport@origene.com EU: info-de@origene.com CN: techsupport@origene.cn

Product datasheet for RC206616L1

Carbonic Anhydrase I (CA1) (NM_001738) Human Tagged Lenti ORF Clone

Product data:

| Product Type: | Expression Plasmids |
|------------------------------|---|
| Product Name: | Carbonic Anhydrase I (CA1) (NM_001738) Human Tagged Lenti ORF Clone |
| Tag: | Myc-DDK |
| Symbol: | Carbonic Anhydrase I |
| Synonyms: | CA-l; CAB; Car1; HEL-S-11 |
| Mammalian Cell Selection: | None |
| Vector: | pLenti-C-Myc-DDK (PS100064) |
| E. coli Selection: | Chloramphenicol (34 ug/mL) |
| ORF Nucleotide Sequence: | The ORF insert of this clone is exactly the same as(RC206616). |
| Restriction Sites: | Sgfl-Mlul |
| Cloning Scheme: | |
| | Cloning sites used for ORF Shuttling: |
| | Sgf I ORF Mlu I GCG ATC GC C ATG// NNN ACG CGT |

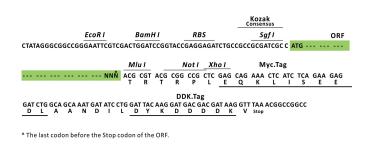

ACCN: ORF Size: NM_001738 783 bp



View online »


This product is to be used for laboratory only. Not for diagnostic or therapeutic use. ©2023 OriGene Technologies, Inc., 9620 Medical Center Drive, Ste 200, Rockville, MD 20850, US

| | nic Anhydrase I (CA1) (NM_001738) Human Tagged Lenti ORF Clone – RC206616L1 |
|------------------------|--|
| OTI Disclaimer: | The molecular sequence of this clone aligns with the gene accession number as a point of reference only. However, individual transcript sequences of the same gene can differ through naturally occurring variations (e.g. polymorphisms), each with its own valid existence. This clone is substantially in agreement with the reference, but a complete review of all prevailing variants is recommended prior to use. <u>More info</u> |
| OTI Annotation: | This clone was engineered to express the complete ORF with an expression tag. Expression varies depending on the nature of the gene. |
| Components: | The ORF clone is ion-exchange column purified and shipped in a 2D barcoded Matrix tube containing 10ug of transfection-ready, dried plasmid DNA (reconstitute with 100 ul of water). |
| Reconstitution Method: | Centrifuge at 5,000xg for 5min. Carefully open the tube and add 100ul of sterile water to dissolve the DNA. Close the tube and incubate for 10 minutes at room temperature. Briefly vortex the tube and then do a quick spin (less than 5000xg) to concentrate the liquid at the bottom. Store the suspended plasmid at -20°C. The DNA is stable for at least one year from date of shipping when stored at -20°C. |
| RefSeq: | <u>NM 001738.1</u> |
| RefSeq Size: | 1319 bp |
| RefSeq ORF: | 786 bp |
| Locus ID: | 759 |
| UniProt ID: | <u>P00915</u> |
| Cytogenetics: | 8q21.2 |
| Domains: | carb_anhydrase |
| Protein Families: | Druggable Genome |
| Protein Pathways: | Nitrogen metabolism |
| MW: | 28.9 kDa |
| Gene Summary: | Carbonic anhydrases (CAs) are a large family of zinc metalloenzymes that catalyze the reversible hydration of carbon dioxide. They participate in a variety of biological processes, including respiration, calcification, acid-base balance, bone resorption, and the formation of aqueous humor, cerebrospinal fluid, saliva and gastric acid. They show extensive diversity in tissue distribution and in their subcellular localization. This CA1 gene is closely linked to the CA2 and CA3 genes on chromosome 8. It encodes a cytosolic protein that is found at the highest level in erythrocytes. Allelic variants of this gene have been described in some populations. Alternative splicing and the use of alternative promoters results in multiple |


This product is to be used for laboratory only. Not for diagnostic or therapeutic use. ©2023 OriGene Technologies, Inc., 9620 Medical Center Drive, Ste 200, Rockville, MD 20850, US

transcript variants. [provided by RefSeq, Nov 2016]

Product images:

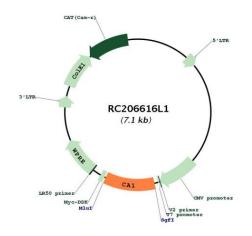

Circular map for RC206616L1

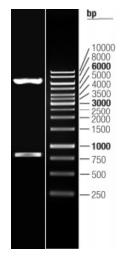

Double digestion of RC206616L1 using Sgfl and Mlul

This product is to be used for laboratory only. Not for diagnostic or therapeutic use. ©2023 OriGene Technologies, Inc., 9620 Medical Center Drive, Ste 200, Rockville, MD 20850, US