


Product datasheet for RC223769L3

ATP6V0A1 (NM_005177) Human Tagged Lenti ORF Clone

Product data:

Product Type: Expression Plasmids Product Name: ATP6V0A1 (NM_005177) Human Tagged Lenti ORF Clone Tag: Myc-DDK Symbol: ATP6V0A1 Synonyms: a1; ATP6N1; ATP6N1A; Stv1; Vph1; VPP1 Mammalian Cell Puromycin Selection: Vector: pLenti-C-Myc-DDK-P2A-Puro (PS100092) E. coli Selection: Chloramphenicol (34 ug/mL) The ORF insert of this clone is exactly the same as(RC223769). **ORF** Nucleotide Sequence: **Restriction Sites:** Sgfl-Mlul **Cloning Scheme:** Cloning sites used for ORF Shuttling: ORF Sqf I Mlu I

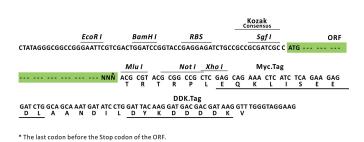

--- GCG ATC GC C ATG --- //--- NNN ACG CGT ---

ACCN: NM_005177 ORF Size: 2493 bp





View online »

This product is to be used for laboratory only. Not for diagnostic or therapeutic use. ©2023 OriGene Technologies, Inc., 9620 Medical Center Drive, Ste 200, Rockville, MD 20850, US

OriGene Technologies, Inc.

9620 Medical Center Drive, Ste 200 Rockville, MD 20850, US Phone: +1-888-267-4436 https://www.origene.com techsupport@origene.com EU: info-de@origene.com CN: techsupport@origene.cn

ORIGENE ATP6VO	DA1 (NM_005177) Human Tagged Lenti ORF Clone – RC223769L3
OTI Disclaimer:	The molecular sequence of this clone aligns with the gene accession number as a point of reference only. However, individual transcript sequences of the same gene can differ through naturally occurring variations (e.g. polymorphisms), each with its own valid existence. This clone is substantially in agreement with the reference, but a complete review of all prevailing variants is recommended prior to use. <u>More info</u>
OTI Annotation:	This clone was engineered to express the complete ORF with an expression tag. Expression varies depending on the nature of the gene.
Components:	The ORF clone is ion-exchange column purified and shipped in a 2D barcoded Matrix tube containing 10ug of transfection-ready, dried plasmid DNA (reconstitute with 100 ul of water).
Reconstitution Method:	 Centrifuge at 5,000xg for 5min. Carefully open the tube and add 100ul of sterile water to dissolve the DNA. Close the tube and incubate for 10 minutes at room temperature. Briefly vortex the tube and then do a quick spin (less than 5000xg) to concentrate the liquid at the bottom. Store the suspended plasmid at -20°C. The DNA is stable for at least one year from date of shipping when stored at -20°C.
RefSeq:	<u>NM 005177.2</u>
RefSeq Size:	4139 bp
RefSeq ORF:	2496 bp
Locus ID:	535
UniProt ID:	<u>Q93050</u>
Cytogenetics:	17q21.2
Domains:	V_ATPase_sub_a
Protein Families:	Transmembrane
Protein Pathways:	Epithelial cell signaling in Helicobacter pylori infection, Lysosome, Metabolic pathways, Oxidative phosphorylation, Vibrio cholerae infection
MW:	95.6 kDa

This product is to be used for laboratory only. Not for diagnostic or therapeutic use. ©2023 OriGene Technologies, Inc., 9620 Medical Center Drive, Ste 200, Rockville, MD 20850, US

Gene Summary: This gene encodes a component of vacuolar ATPase (V-ATPase), a multisubunit enzyme that mediates acidification of eukaryotic intracellular organelles. V-ATPase dependent organelle acidification is necessary for such intracellular processes as protein sorting, zymogen activation, receptor-mediated endocytosis, and synaptic vesicle proton gradient generation. V-ATPase is composed of a cytosolic V1 domain and a transmembrane V0 domain. The V1 domain consists of three A and three B subunits, two G subunits plus the C, D, E, F, and H subunits. The V1 domain contains the ATP catalytic site. The V0 domain consists of five different subunits: a, c, c', c'', and d. Additional isoforms of many of the V1 and V0 subunit proteins are encoded by multiple genes or alternatively spliced transcript variants. This gene encodes one of three A subunit proteins and the encoded protein is associated with clathrin-coated vesicles. Three transcript variants encoding different isoforms have been found for this gene. [provided by RefSeq, Jul 2008]

This product is to be used for laboratory only. Not for diagnostic or therapeutic use. ©2023 OriGene Technologies, Inc., 9620 Medical Center Drive, Ste 200, Rockville, MD 20850, US