

#### OriGene Technologies, Inc.


9620 Medical Center Drive, Ste 200 Rockville, MD 20850, US Phone: +1-888-267-4436 https://www.origene.com techsupport@origene.com EU: info-de@origene.com CN: techsupport@origene.cn

# Product datasheet for RC226702L4

#### PKR (EIF2AK2) (NM\_001135652) Human Tagged Lenti ORF Clone

### **Product data:**

| Product Type:                | Expression Plasmids   |
|------------------------------|---|
| Product Name:                | PKR (EIF2AK2) (NM_001135652) Human Tagged Lenti ORF Clone                                   |
| Tag:                         | mGFP  |
| Symbol:                      | PKR   |
| Synonyms:                    | EIF2AK1; LEUDEN; PKR; PPP1R83; PRKR   |
| Mammalian Cell<br>Selection: | Puromycin   |
| Vector:                      | pLenti-C-mGFP-P2A-Puro (PS100093)   |
| E. coli Selection:           | Chloramphenicol (34 ug/mL)  |
| ORF Nucleotide<br>Sequence:  | The ORF insert of this clone is exactly the same as(RC226702).                              |
| <b>Restriction Sites:</b>    | Sgfl-Mlul   |
| Cloning Scheme:              |   |
|                              | Cloning sites used for ORF Shuttling:   |
|                              | Sgf I         ORF         Mlu I            GCG ATC GC C         ATG / / NNN         ACG CGT |

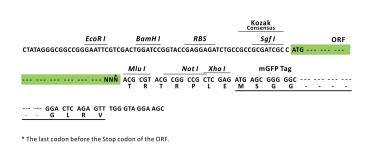

ACCN: ORF Size: NM\_001135652 1530 bp




View online »

This product is to be used for laboratory only. Not for diagnostic or therapeutic use. ©2023 OriGene Technologies, Inc., 9620 Medical Center Drive, Ste 200, Rockville, MD 20850, US

| PKR (EIF2AK2) (NM_001135652) Human Tagged Lenti ORF Clone – RC226702L4 |  |
|--|--|
| OTI Disclaimer:  | The molecular sequence of this clone aligns with the gene accession number as a point of reference only. However, individual transcript sequences of the same gene can differ through naturally occurring variations (e.g. polymorphisms), each with its own valid existence. This clone is substantially in agreement with the reference, but a complete review of all prevailing variants is recommended prior to use. <u>More info</u>                                  |
| OTI Annotation:  | This clone was engineered to express the complete ORF with an expression tag. Expression varies depending on the nature of the gene.   |
| Components:  | The ORF clone is ion-exchange column purified and shipped in a 2D barcoded Matrix tube containing 10ug of transfection-ready, dried plasmid DNA (reconstitute with 100 ul of water).   |
| Reconstitution Method:   | <ol> <li>Centrifuge at 5,000xg for 5min.</li> <li>Carefully open the tube and add 100ul of sterile water to dissolve the DNA.</li> <li>Close the tube and incubate for 10 minutes at room temperature.</li> <li>Briefly vortex the tube and then do a quick spin (less than 5000xg) to concentrate the liquid at the bottom.</li> <li>Store the suspended plasmid at -20°C. The DNA is stable for at least one year from date of shipping when stored at -20°C.</li> </ol> |
| RefSeq:  | <u>NM 001135652.1, NP 001129124.1</u>  |
| RefSeq ORF:  | 1533 bp  |
| Locus ID:  | 5610   |
| UniProt ID:  | <u>P19525</u>  |
| Cytogenetics:  | 2p22.2   |
| Protein Families:  | Druggable Genome, Protein Kinase, Transcription Factors  |
| MW:  | 57.2 kDa   |
| Gene Summary:  | The protein encoded by this gene is a serine/threonine protein kinase that is activated by autophosphorylation after binding to dsRNA. The activated form of the encoded protein can phosphorylate translation initiation factor EIF2S1, which in turn inhibits protein synthesis. This protein is also activated by manganese ions and heparin. Three transcript variants encoding two different isoforms have been found for this gene. [provided by RefSeq, Oct 2011]   |

This product is to be used for laboratory only. Not for diagnostic or therapeutic use. ©2023 OriGene Technologies, Inc., 9620 Medical Center Drive, Ste 200, Rockville, MD 20850, US

## **Product images:**

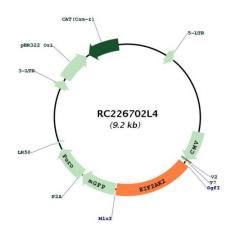

Circular map for RC226702L4

This product is to be used for laboratory only. Not for diagnostic or therapeutic use. ©2023 OriGene Technologies, Inc., 9620 Medical Center Drive, Ste 200, Rockville, MD 20850, US