

OriGene Technologies, Inc.

9620 Medical Center Drive, Ste 200 Rockville, MD 20850, US Phone: +1-888-267-4436 https://www.origene.com techsupport@origene.com EU: info-de@origene.com CN: techsupport@origene.cn

Product datasheet for UM570025

PSMA (FOLH1) Mouse Monoclonal Antibody [Clone ID: UMAB25]

Product data:

Product Type:	Primary Antibodies
Clone Name:	UMAB25
Applications:	10k-ChIP, IF, IHC, WB
Recommended Dilution:	WB 1:500~2000, IF 1:100, IF 1:100
Reactivity:	Human, Dog, Rat, Monkey, Mouse
Host:	Mouse
lsotype:	lgG2b
Clonality:	Monoclonal
Immunogen:	Full length human recombinant protein of human FOLH1(NP_004467) produced in HEK293T cell.
Formulation:	PBS (pH 7.3) containing 1% BSA, 50% glycerol and 0.02% sodium azide.
Concentration:	0.5~1.0 mg/ml (Lot Dependent)
Purification:	Purified from mouse ascites fluids or tissue culture supernatant by affinity chromatography (protein A/G)
Conjugation:	Unconjugated
Storage:	Store at -20°C as received.
Stability:	Stable for 12 months from date of receipt.
Predicted Protein Size:	84.2 kDa
Gene Name:	folate hydrolase 1
Database Link:	<u>NP 004467</u> Entrez Gene 53320 MouseEntrez Gene 85309 RatEntrez Gene 476775 DogEntrez Gene 707714 MonkeyEntrez Gene 2346 Human Q04609




SMA (FOLH1) Mouse Monoclonal Antibody [Clone ID: UMAB25] – UM570025 PSMA (FOLH1) Mouse Monoclonal Antibody [Clone ID: UMAB25] – UM570025

Background: This gene encodes a type II transmembrane glycoprotein belonging to the M28 peptidase family. The protein acts as a glutamate carboxypeptidase on different alternative substrates, including the nutrient folate and the neuropeptide N-acetyl-l-aspartyl-l-glutamate and is expressed in a number of tissues such as prostate, central and peripheral nervous system and kidney. A mutation in this gene may be associated with impaired intestinal absorption of dietary folates, resulting in low blood folate levels and consequent hyperhomocysteinemia. Expression of this protein in the brain may be involved in a number of pathological conditions associated with glutamate excitotoxicity. In the prostate the protein is upregulated in cancerous cells and is used as an effective diagnostic and prognostic indicator of prostate cancer. This gene likely arose from a duplication event of a nearby chromosomal region. Alternative splicing gives rise to multiple transcript variants encoding several different isoforms. [provided by RefSeq]


Synonyms: FGCP; FOLH; GCP2; GCPII; mGCP; NAALAD1; NAALAdase; PSM; PSMA

Protein Families: Druggable Genome, Protease, Transmembrane

Product images:

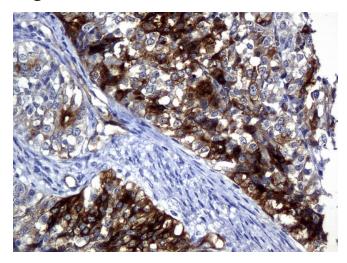


Immunohistochemical staining of paraffinembedded Adenocarcinoma of endometrium tissue using anti-FOLH1mouse monoclonal antibody. (Clone UMAB25, dilution 1:100; heatinduced epitope retrieval by 10mM citric buffer, pH6.0, 120°C for 3min)

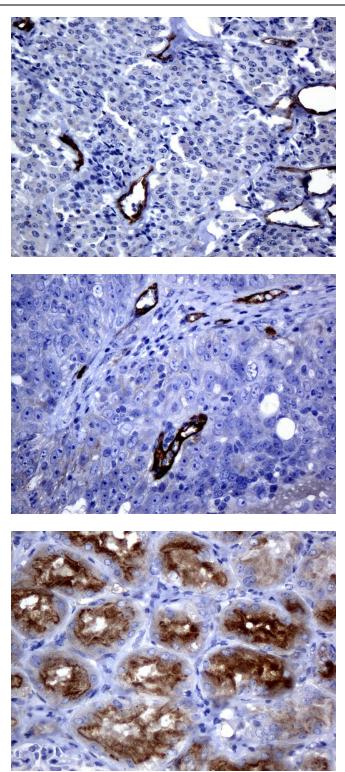

Immunohistochemical staining of paraffinembedded Carcinoma of thyroid tissue using anti-FOLH1mouse monoclonal antibody. (Clone UMAB25, dilution 1:100; heat-induced epitope retrieval by 10mM citric buffer, pH6.0, 120°C for 3min)

Immunohistochemical staining of paraffinembedded Adenocarcinoma of ovary tissue using anti-FOLH1mouse monoclonal antibody. (Clone UMAB25, dilution 1:100; heat-induced epitope retrieval by 10mM citric buffer, pH6.0, 120°C for 3min)


Immunohistochemical staining of paraffinembedded Kidney tissue using anti-FOLH1mouse monoclonal antibody. (Clone UMAB25, dilution 1:100; heat-induced epitope retrieval by 10mM citric buffer, pH6.0, 120°C for 3min)

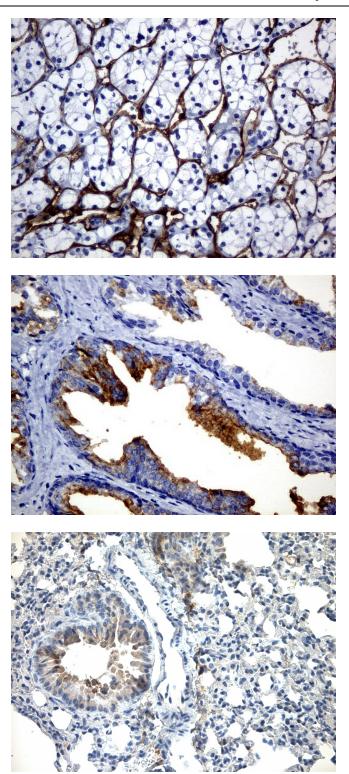


Immunohistochemical staining of paraffinembedded Carcinoma of kidney tissue using anti-FOLH1mouse monoclonal antibody. (Clone UMAB25, dilution 1:100; heat-induced epitope retrieval by 10mM citric buffer, pH6.0, 120°C for 3min)

Immunohistochemical staining of paraffinembedded Carcinoma of prostate tissue using anti-FOLH1mouse monoclonal antibody. (Clone UMAB25, dilution 1:100; heat-induced epitope retrieval by 10mM citric buffer, pH6.0, 120°C for 3min)


Immunohistochemical staining of paraffinembedded mouse lung tissue using anti-FOLH1 (PSMA) clone UMAB25 mouse monoclonal antibody. HIER TEE buffer pH9 ([B21-100]) at 110C for 10 min, [UM500025] (1:100). Detection was done with Klear Mouse (C/N [D52-18]) DAB Kit.

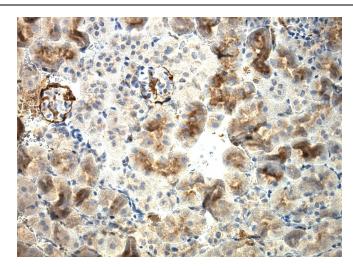


Immunohistochemical staining of paraffinembedded mouse kidney tissue using anti-FOLH1 (PSMA) clone UMAB25 mouse monoclonal antibody. HIER TEE buffer pH9 ([B21-100]) at 110C for 10 min, [UM500025] (1:100). Detection was done with Klear Mouse (C/N [D52-18]) DAB Kit.

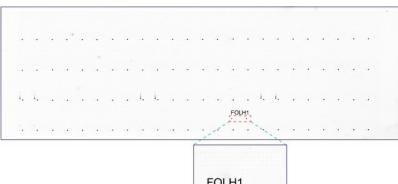


OriGene overexpression protein microarray chip was immunostained with UltraMAB anti-FOLH1 mouse monoclonal antibody (Clone UMAB25). The positive reactive proteins are highlighted with two red arrows in the enlarged subarray. All the positive controls spotted in this subarray are also labeled for clarification. These data show that UltraMAB anti-FOLH1 (Clone UMAB25) very specifically recognizes FOLH1 antigen on OriGene protein microarray chip.

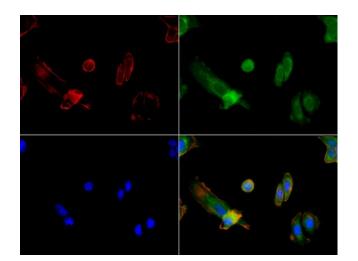


Immunofluorescent staining of PC-12 cells using FOLH1 mouse monoclonal antibody ([UM500025], green). Actin filaments were labeled with TRITC-phalloidin (red), and nuclear with DAPI (blue). The three-color overlay image is located at the bottom-right corner.

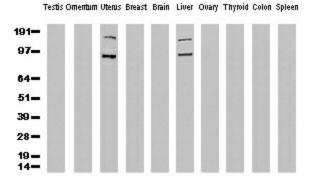


Western Blot analysis of 10 different human tissue lysates (10ug) by using anti-FOLH1 monoclonal antibody (clone UMAB25, 1:500)

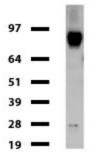


Western blot of cell lysates (35ug) from MDCK. Diluation: 1:500

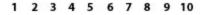

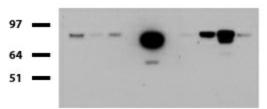

Western blot of human tissue lysates (15ug) from 10 different tissues (1: Testis, 2: Omentum, 3: Uterus, 4: Breast, 5: Brain, 6: Thyroid, 7: Colon, 8: Spleen 9: Liver, 10: Ovary). Diluation: 1:500.

This product is to be used for laboratory only. Not for diagnostic or therapeutic use. ©2024 OriGene Technologies, Inc., 9620 Medical Center Drive, Ste 200, Rockville, MD 20850, US

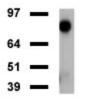

Western blot of mouse tissue lysates (20ug) from Brain. Diluation: 1:500.

This product is to be used for laboratory only. Not for diagnostic or therapeutic use. ©2024 OriGene Technologies, Inc., 9620 Medical Center Drive, Ste 200, Rockville, MD 20850, US